



MASSAGE GUN MB5 MINI

USER MANUAL [ENGLISH](#)



Model: P24MB5GR, P24MB5WH, P24MB5PU



Powerful
Brushless Motor



Active
Lower noise



Light weight
and Compact



C-Type
Rechargeable

PRODUCT OVERVIEW



PRODUCT LIST

- 1x | The device
- 1x | Type-C charging cable
- 1x | The instructions
- 1x | Bullet Head
- 1x | Ball Head
- 1x | Flat Head
- 1x | Fork Head

PRODUCT PARAMETERS

Model: MB5 MINI massage Gun

Rated voltage: 7.4v

Rated current: 2A


Rated power: 25W

Input parameter: 5V  2A

Maximum speed: 3200r/min

Battery capacity: 800mAh


PRODUCT USE

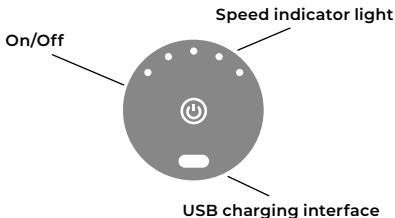

On/Off

Hold the  key for three seconds to turn on

Hold the  key for three seconds to turn off


Click the  button to adjust the speed

The maximum speed is five gear


The starting gear is the lowest gear

Note: in order to protect the bones and muscles, avoid long-term fixed-point massage, resulting in damage.

USING THE MESSAGE HEADS

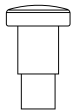


Apply to all
parts
of the body

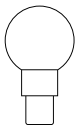

Apply to large
muscle
groups




Apply to all
joints
of the body

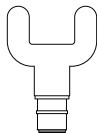


Apply
to the
spine

CHANGE MESSAGE HEAD

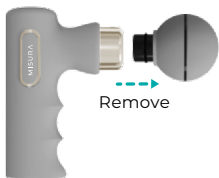

CHARGING

When the power is less than 5% the indicator will flash.
At this time, the machine should be charging.

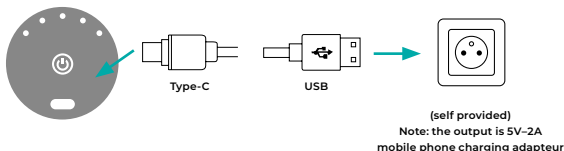

- Please be sure to fully charge before first use.
- The indicator is flashing when charging, and bright when charging finish.
- The charging time is about 4 hours.
- To ensure the safety of use, do not use the massage gun when charging.
- To protect the battery life, do not wait for the battery to run out of charge.

INDICATION

- Pain and spasm due to muscle contusion, sprain and strain.
- To help the swelling flow.
- Loose and thickened connective tissue and fascia.
- Reduce the accumulation of lactic acid in muscle.
- Increase joint mobility.
- Eliminate muscle fatigue.

SAFETY INSTRUCTIONS

To reduce the risk of electric shock, fire or personal injury, this product must be used in accordance with the following instructions.

- Adult use only.
- Use only on a dry clean body surface or over clothing by gently pressing and moving over the body for approximately 60 seconds in each area.
- Use the massager only on soft body tissue as needed without causing pain or discomfort.
- Do not use it on the head or on hard or bony areas of the body.
- Use the heads that give the best desired effect on the specific muscle area.
- Bruising may occur regardless of the control setting or pressure applied. Check frequently treated areas and stop the device immediately at the first sign of pain or discomfort.
- Keep fingers, hair, or other body parts well away from the shaft and back of the applicator, as pinching may occur.
- Do not insert any objects into the ventilation holes of the massager.
- Do not immerse the device in water and do not allow water to enter the massage device vents.
- Protect the device from falling or misuse.
- Carefully check the device and battery before each use. Do not tamper with or alter the individual components of the massager in any way.
- Never leave the massager unattended when in use or charging.

WARNINGS

- Do not immerse the device in water.
- Protect from liquids and heat sources.
- Protect the ventilation openings from dust and small objects.
- Do not remove screws or attempt to repair the device yourself. Any tampering with the product design will void the warranty on this product.
- Do not use continuously for more than one hour.
- Allow the device to rest for at least 30 minutes before using it again.
- Disconnect the charger from the device battery after charging or before use.
- Do not injure yourself. Consult your physician if you are uncertain about the suitability of using the device.
- Use the product only in accordance with the instructions for use provided in the product manual.

NOTICE

Do not use massage equipment or any vibrating device without prior approval from your doctor, if any of the following applies to you:

Pregnancy, diabetes with complications such as neuropathy or retinal damage, heart weakness, recent surgery, epilepsy and migraines, herniated discs, spondylolisthesis or spondylosis, recent joint replacements or IUDs, metal pins or plates, or any concerns about your physical health. Fragile individuals or children should be accompanied by an adult when using any vibration equipment. These contraindications do not mean that you cannot use vibration or exercise equipment, but consult a physician

first. Ongoing research into the effects of vibration exercise on specific medical disorders will likely lead to a shortening of the list of contraindications listed above. Practical experience has shown that integrating vibration exercises into a treatment plan is recommended in several cases. The exercises must be carried out on the advice and with the assistance of a physician or other specialist or physiotherapist.

WARRANTY

The warranty runs from the time of purchase. If the product fails due to defects in materials or design process within 24 months, the Company will repair or replace the defective parts or replace new products free of charge, except for the following:

- Misuse or damage to the equipment caused by shipping.
- Unauthorized disassembly or repair of this equipment.
- Failure to follow the instructions for use.
- Product is damaged due to unusual storage or maintenance environment of the customer.
- If proof of purchase is not provided, the company reserves the right to refuse warranty.

LEGAL DISCLAIMER

If you suffer from any of the conditions listed in the Cautions section, always consult your physician when using MISURA MB5 MINI. The MISURA Store s.r.o. is not a licensed health care provider and represents that it has no expertise in diagnosing, investigating or treating medical conditions of any kind or in determining the effect of any particular exercise or therapy on any medical condition. There is a possibility of physical injury when participating in any exercise or therapy. If you engage in such exercise or therapy, you agree to do so at your own risk, voluntarily participate in these activities, and thereby assume all risk of injury to yourself. By doing so, you also agree to release any and all claims, known or unknown, against MISURA Store s.r.o. arising out of negligence or misuse of the product.

You must be 18 years of age or older. If you experience severe pain or discomfort, stop using the device immediately. Any exercise program, including this one, may result in injury. Consult a physician before starting this or any exercise or therapy program. The manufacturers, distributors, suppliers and dealers disclaim any liability for any losses in connection with the information, exercise demonstrations or instructions contained herein.

Copyright© 2024. MISURA Store s.r.o., 24 month limited warranty. Validity of Terms and Conditions.



Thank you for your purchase.
If necessary, contact us at:

Producer / Importer

MISURA Store s.r.o.,
ID: CZ264 03 714
www.misura.shop
info@misura.shop



We are not liable for any typographical errors.